Restek ProFLOW 6000 Electronic Flowmeter

Version 5.6
(Cat.# 22656)

Table of Contents

- 1.0 Introduction 3
- 2.0 Specifications..... 4
- 3.0 Installing the batteries 4
- 4.0 Battery power consumption 5
 - 4.1 Battery lifetime..... 5
 - 4.2 Battery charge indicator 5
- 5.0 Operating instructions 5
- 6.0 Interpreting results 6
 - 6.1 Flow range display..... 6
- 7.0 Data collection on the PC 7
- 8.0 ProFLOW 6000 flowmeter menu structure..... 9
 - 8.1 Unit power up/power down..... 9
 - 8.2 USB activation..... 9
 - 8.3 Adjust LCD image backlight..... 10
 - 8.4 Show battery charge indicator..... 10
 - 8.5 Firmware version information 11
 - 8.6 Adjust auto shutoff duration..... 11
 - 8.7 Display calibration date 12
 - 8.8 Display unit serial number 12
- 9.0 Troubleshooting 12
- 10.0 Product back label legend 13
- 11.0 Volumetric vs. mass flow measurements 14
- 12.0 Bubble flowmeters 14
- 13.0 Product case specifications 14
- 14.0 Calibration and service 14
 - Declaration of Conformity 31

Restek ProFLOW 6000 Electronic Flowmeter

1.0 Introduction


Restek's ProFLOW 6000 flowmeter is specifically designed for use with gas chromatography (GC) systems. The probe is applied directly to the gas flow stream, and the measured flow rate is presented on the LCD screen. Units of flow are measured in mL/min.

This unit provides continuous real-time measurements of gas streams ranging from 0.50 mL/min to 500 mL/min. Because the technology uses volumetric flow measurement, the unit is compatible with all laboratory gases.

CAUTION: Do NOT exceed maximum operating flow rates. Recalibration may be required if the unit has been subjected to extreme flow rates.

Always use appropriate laboratory safety practices when operating this device. Wear laboratory safety goggles when operating this unit.

Figure 1: Product diagrams

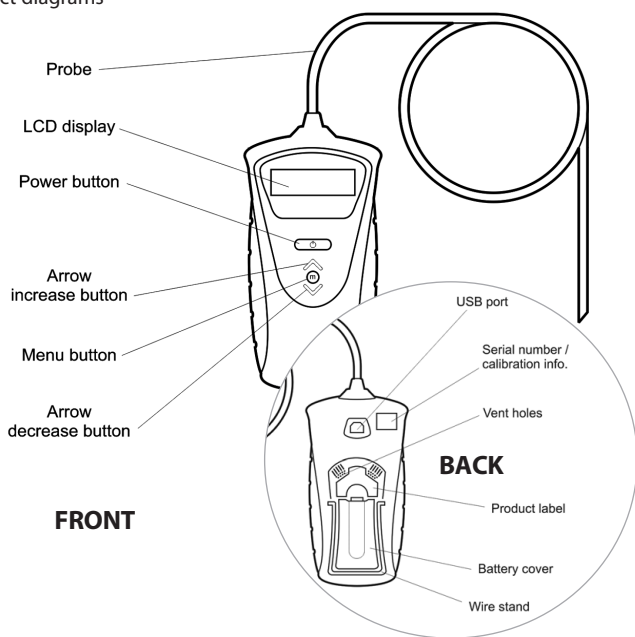NOTE: There are no serviceable parts in this unit. Opening the device—other than to change the batteries—or tampering with the internal parts will void the factory warranty.

NOTE: To ensure accurate measurements and effective clearance of the flow gas from the unit, DO NOT obstruct the vent holes on the back of the unit.

NOTE: Modifying the length of tubing can cause a shift in calibration.

2.0 Specifications

Table I:

Type of measurement	Volumetric flow
Accuracy of measurements	Accuracy of ± 2.00% of flow reading or ± 0.200 mL/min, whichever is greater.
Power requirements	2 AA alkaline batteries 1.5 V DC each/3 V DC 200 ma
Operating flow range	0.50 to 500 mL/min
Operating temperature range	32 °-120 °F (0 °- 48 °C)
Available communication	USB type B data port
Warranty	One year (excludes recalibration)
Calibration	NIST traceable Yearly recalibration is recommended
Certification/Compliance	CE, Ex, Canadian ICES-003, WEEE, RoHS 2, China RoHS 2, UKCA (see section 10.0)

3.0 Installing the batteries

This unit uses two AA alkaline batteries.

To install batteries, extend the wire stand. Open the cover. Insert the batteries with the polarity (⊕ and ⊖) correctly aligned. Close the cover. (Figures 2 and 3.)

Precautions for battery replacement:

- Load the new batteries with their polarity (⊕ and ⊖) aligned correctly.
- **Do not use rechargeable batteries.**


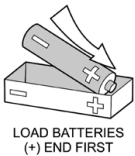
Figure 2: Extend the wire stand before opening battery door.

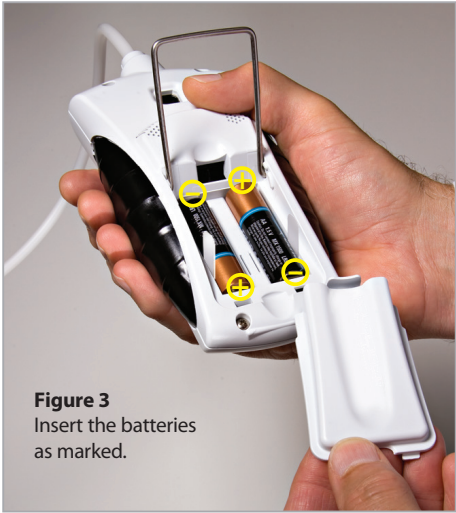4.0 Battery power consumption

4.1 Battery lifetime

The battery lifetime is dependent on the number of options the user has enabled.

The unit is shipped with the most energy-demanding options disabled (Table II).

The power saving functions can be changed.

- ➔ See Section 8.2: USB Activation
- ➔ See Section 8.3: Adjust LCD Image Backlight
- ➔ See Section 8.6: Adjust Auto Shutoff Duration

Table II: Default settings for the ProFLOW 6000 flowmeter

Auto shutoff duration	6 minutes
LCD backlight	0 (off)
USB port	disabled

NOTE: Store your ProFLOW 6000 flowmeter in its protective storage case following use. Do not store this manual or any other items on top of the ProFLOW 6000 flowmeter or the unit may turn on when the storage case lid is closed.

4.2 Battery charge indicator

The unit includes a battery charge indicator. Replace batteries as needed.


- ➔ See Section 8.4: Show Battery Charge Indicator

5.0 Operating instructions

 **CAUTION: Do not exceed maximum operating flow rates. Recalibration may be required if the unit has been subjected to extreme flow rates.**

Connect the white probe end tip to the output of the gas flow line to be measured. Be sure the probe tip connection is completely sealed around the flow source outlet and is free of leaks (Figure 4). (White tubing ID 0.125"; OD 0.250")

Press and hold the  (Power) button until the unit responds with a regular clicking sound. The ProFLOW will immediately begin to provide flow measurements (Figure 5). Wait for the measured values to stabilize. It takes a few seconds for the unit to reach a steady state with the gas flow line.

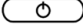To power down the unit, press and hold the  (Power) button until the unit stops clicking.

The unit is equipped with a timed auto shutoff option. (Default: 6 minutes.)

- ➔ See Section 8.6: Adjust Auto Shutoff Duration

Figure 4: Probe connected to a GC gas outlet.

Figure 5: LED displays the measured flow value.

6.0 Interpreting results

The unit has an operating range of 0.50 mL/min to 500 mL/min (Figure 7). If the flow is less than 0.50 mL/min., the display will read “under range.” If the flow exceeds 515 mL/min, the display will read “over range.” Excessively high flow rates may damage this unit.

NOTE: units of mL/min are equivalent to ccm.

6.1 Flow range display

The unit automatically adjusts the resolution of the display depending on the flow range being measured. Table III shows the resolution of the flow ranges.


Table III: Display resolution vs. flow range.

Flow range	Display resolution (mL/min)
0.50–9.99	0.01
10.0–99.9	0.1
100–500	1

Figure 6: Example flow value.

Figure 7: Description of flow ranges.

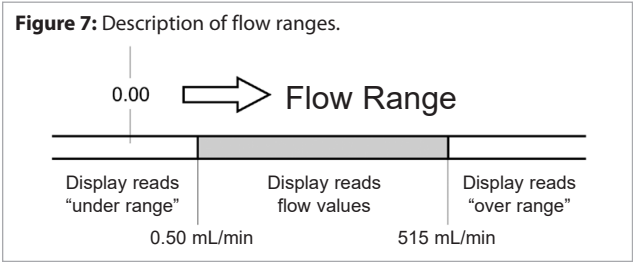7.0 Data collection on the PC

WARNING: ONLY connect USB cable to USB port while unit is OFF.

The ProFLOW 6000 flowmeter provides you with a data stream of real-time flow values via the USB port (Figure 1, p. 3). In order to use this feature, you must first install the appropriate FTDI Virtual Com Port (VCP) Driver available at:

<http://www.ftdichip.com/Drivers/VCP.htm>


The VCP driver will cause the ProFLOW 6000 flowmeter to appear as a standard RS-232 port. This will work on any operating system for which there is an FTDI VCP driver. After installing the driver, connecting the device, and determining which port it creates, you can access the data stream through any programmatic means, or by using any serial terminal software.

For Windows systems:

To determine which port the ProFLOW 6000 flowmeter is using, go to the Control Panel and open System. Go to the Hardware tab and click the Device Manager button. Expand the Ports (COM & LPT) entry. Make sure the VCP driver is installed and then connect a powered ProFLOW 6000 flowmeter to the USB port. You will see the new COM port appear. Open your serial terminal.

➔ See Section 8.2: USB Activation

Figure 8: Screen capture of data collection.

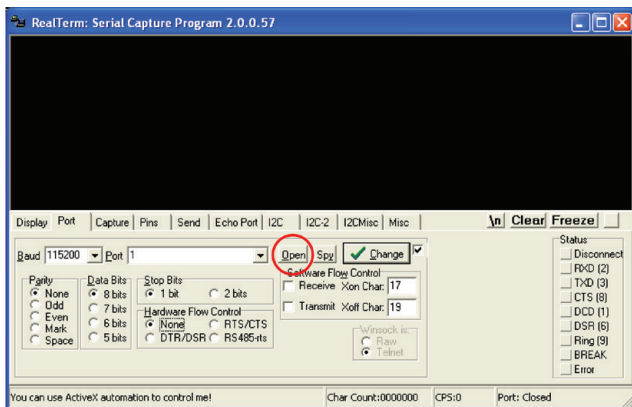If you do not already have serial terminal software, free, open-source options are available online (i.e., RealTerm, etc.). You can download RealTerm* software from:

<http://realterm.sourceforge.net/>

After installation, click the Port tab and set the following:


- Baud: 115200
- Port: The appropriate VCP for your ProFLOW 6000 flowmeter
- Parity: None
- Data Bits: 8 bits
- Stop Bits: 1 bit
- HardwareFlow Control: None

Go to the USB menu entry in the ProFLOW 6000 flowmeter and turn the transmission on. Finally, click **Open** on the Port tab in RealTerm and you will see the serial data stream begin in the terminal window (Figure 8).

If you would like to log the flow data, this can be done by clicking on the Capture tab. Set File to the name and location of the log file that you would like to save and click either the **Start: Overwrite** or **Start: Append** buttons appropriately (Figure 9).

*This software is not supplied or supported by Restek. User assumes all responsibility for the downloading and use of the program.

Figure 9: Screen capture of data collection.

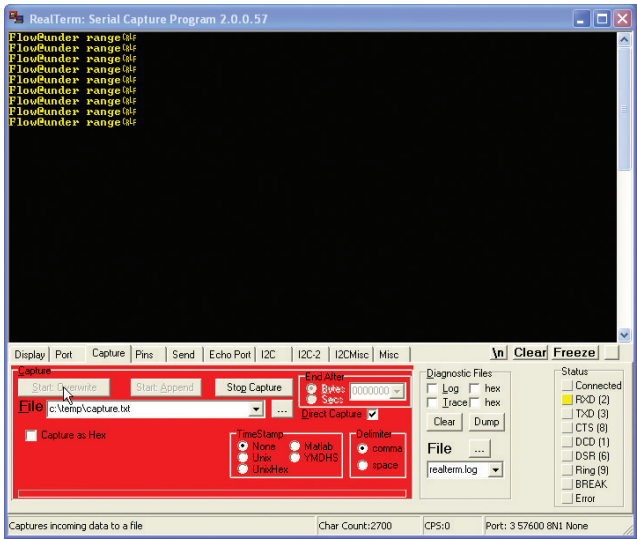8.0 ProFLOW 6000 flowmeter menu structure

8.1 Unit power up/power down

Press the  (Power) button:

The LCD screen will display the device intro screen:

Followed by calibrated date:

01/10/2012

Followed by measured flow data:

50.5 ml/min

To power off, press and hold the  (Power) button.

Other messages encountered at power up

under range

If the unit is hooked up to a flow stream with a flow rate less than 0.50 mL/min, the unit will report an “under range” status for the flow. This message will appear until the flow rate exceeds 0.50 mL/min.

8.2 USB activation

To enable the USB, press the  (Menu) button.

Use the   (arrow) keys to select the USB menu:

►USB

Press the  (menu) button again to enter the value select screen.

Use the   (arrow) keys to toggle between USB “on” and “off.”

To return to the main menu screen, press the  (menu) button.


To exit and return to measuring flow, press the  (menu) button again.

➔ See Section 7.0: Data Collection on the PC

8.3 Adjust LCD image backlight

Press the  (Menu) button.

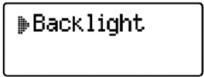Use the   (arrow) keys to select the Backlight menu:

Press the  (Menu) button again to enter the value select screen.

Use the   (arrow) keys to select the backlight value.

Backlight Values: 0 (off) and 5 (maximum).


To return to the main menu screen, press the  (Menu) button.

To exit and return to measuring flow, press the  (Menu) button again.

8.4 Show battery charge indicator

Press the  (Menu) button.

Use the   (arrow) keys to select the Batteries menu:

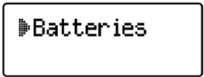Press the  (Menu) button again.

The battery life is displayed.


To return to the main menu screen, press the  (Menu) button.

To exit and return to measuring flow, press the  (Menu) button again.

8.5 Firmware version information


Press the  (Menu) button.

Use the   (arrow) keys to select the Firmware menu:

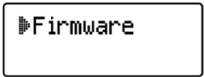Press the  (Menu) button again.

The most recent version of firmware is displayed.

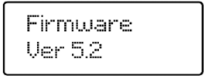
To return to the main menu screen, press the  (Menu) button.

To exit and return to measuring flow, press the  (Menu) button again.

8.6 Adjust auto shutoff duration

To conserve battery life, the unit automatically turns off after 6 minutes. To customize the auto shutoff setting, press the  (Menu) button.

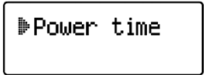Use the   (arrow) keys to select the Power Time menu:

Press the  (Menu) button again to enter the value select screen.

Use the   (arrow) keys to select the auto shutoff setting.

Values: 1–59 minutes or “constant on” (max.)

To return to the main menu screen, press the  (Menu) button.

To exit and return to measuring flow, press the  (Menu) button again.


8.7 Calibration Date

To display the calibration date, press the  (Menu) button.

Use the   (arrow) keys to select the Cal Date menu:

Press the  (Menu) button again to display the calibration date.

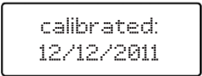
To return to the main menu screen, press the  (Menu) button.

To exit and return to measuring flow, press the  (Menu) button again.


8.8 Serial Number

To display the flowmeter serial number, press the  (Menu) button.

Use the   (arrow) keys to select the Serial Num menu:

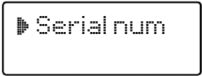Press the  (menu) button again to display the serial number.

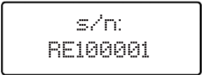To return to the main menu screen, press the  (Menu) button.

To exit and return to measuring flow, press the  (Menu) button again.


9.0 Troubleshooting

Problem	Possible Cause(s)	Suggested Solution(s)
Multiple readings are not giving reproducible results	<ul style="list-style-type: none">• Unit is out of calibration• Value is being compared to a bubble flowmeter	<ul style="list-style-type: none">• Return the unit to Restek for recalibration*• See Section 12.0 for a discussion of the weaknesses of bubble flowmeters
Unit does not power up	<ul style="list-style-type: none">• Dead batteries	<ul style="list-style-type: none">• Replace with 2 new AA alkaline batteries
Flow value display is erratic/jumpy	<ul style="list-style-type: none">• The ProFLOW 6000 flowmeter is very sensitive to small changes in flow	<ul style="list-style-type: none">• Allow more time for flow to stabilize

*Contact returns@restek.com for return instructions for servicing a damaged unit. Additional charges may apply if the warranty has expired or the unit is damaged due to misuse.
Call Restek Customer Service at 1-800-356-1688 or 1-814-353-1300, ext. 3 (or your Restek representative) if you have any questions about this product or any other Restek product.

10.0 Product back label legend

#	Description
1	Serial number
2	Company name
3	This unit conforms to EU/EMC Directive 2004/108/EC; standards to which conformity is declared include 61326:1997 w/A3 Class A.
4	This Class A digital apparatus complies with Canadian ICES-003.
5	This complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
6	Country of origin
7	Product catalog number
8	Company address
9	This unit is WEEE compliant.
10	This unit is UKCA compliant.
11	This unit is RoHS 2 and China RoHS 2 compliant. 
12	
Ex nL	EN60079-0: 2006; Electrical apparatus for explosive gas atmospheres—Part 0: General Requirements. EN60079-15: 2005; Electrical apparatus for explosive gas atmospheres—Part 15: Construction, test, and marking of type of protection “nL” energy limited apparatus.
IIC	Group II applies to areas above ground environments. Gas Group IIC relates to hydrogen and related gas types.
3G	Category 3 relating to gas analysis; normal safety measure. Sufficient safety during normal operation. Normal operation described as measuring flows of flammable or explosive gases in a nonflammable environment.
T6	During testing, neither internal nor external elements exceed 85 °C.
X	Additional information: Operating range: 32 °F ≤ Tamb ≤ 120 °F 0 °C ≤ Tamb ≤ 48 °C Not intended for outdoor use or wet locations.
2010	Year of product design release.

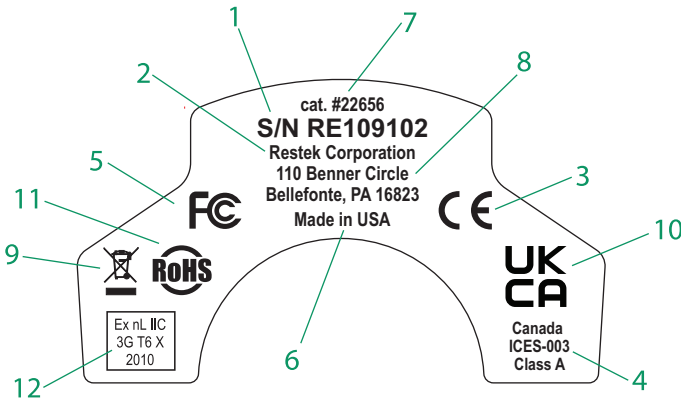11.0 Volumetric vs. mass flow measurements


The Restek ProFLOW 6000 flowmeter is a volumetric flow measurement device. Volumetric flow is the measurement of the volume of gas through a conveyance per quantity of time. Standard units of measure for this parameter are given in mL/min. The advantage of measuring volumetric flow is its independence to the composition of the flow gas. It is not necessary to correct the flow values based on the gas composition, as is required for mass flow devices.

Mass flow measures the weight of the gas flowing through the instrument per quantity of time. Mass flow units of measure are commonly g/sec.

12.0 Bubble flowmeter measurements

If you employ bubble flowmeters in your laboratory, you may find that they give slightly different flow rate values than the Restek ProFLOW 6000 flowmeter. This error is due to technology limitations inherent in the bubble flowmeter device, error from variances in air humidity within the bubble chamber, and its direct contribution to the measured flow rate. In the event a bubble flowmeter is used to measure flow gas where the gas is at elevated temperatures, the error due to humidity contributions can be extreme. For the most accurate measurement of laboratory gas flow rates, we recommend using the Restek ProFLOW 6000 flowmeter over bubble flowmeters.

13.0 Product case specifications

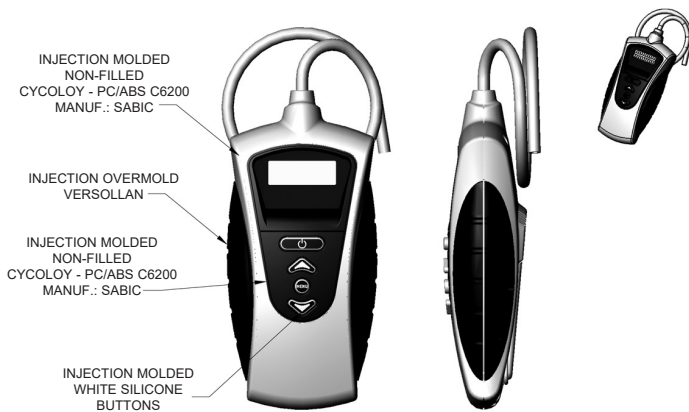14.0 Calibration and service

The Restek ProFLOW 6000 flowmeter comes factory calibrated and carries a one-year warranty (excluding recalibration) from time of purchase. All units are calibrated to NIST traceable standards.

Recommended schedule for recalibration is once every year from time of purchase. Customers will need to return the unit to Restek for recalibration. At that time, preventative maintenance services can also be performed. A fee will be charged for recalibration and servicing of the unit. Prolonged failure to recalibrate the instrument may result in increased error.

Please have the serial number available when calling Restek with any concerns you may have.

Additional charges may apply if the warranty is expired or the damage is due to misuse. This manual is also available in electronic format at www.restek.com.

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DECLARATION OF CONFORMITY

Application of
Council Directive(s):

EU/EMC Directive 2004/108/EC
ATEX Directive 94/9/EC

Standards to which
Conformity is Declared:

61326:1997 w/A3 Class A
ENG1010-1:2001
EN60079-0:2006
EN60079-15:2005
UKCA

Declarer's Name:
Declarer's Address:

Restek Corporation
110 Benner Circle
Bellefonte, PA 16823

Type of Equipment:

ProFLOW 6000 Series
Flow Meter

Catalog number(s)

22656

Person responsible for
making this Declaration:
Position Title:

Gerald Johnston
Product Development Engineer

Date:

06/18/2021

Signature:

Questions about this or any other Restek product? Contact us or your local Restek representative (www.restek.com/contact-us).

Restek patents and trademarks are the property of Restek Corporation. (See www.restek.com/Patents-Trademarks for full list.) Other trademarks in Restek literature or on its website are the property of their respective owners. Restek registered trademarks are registered in the U.S. and may also be registered in other countries.

© 2021 Restek Corporation. All rights reserved. Printed in the U.S.A.

www.restek.com

#204-04-002 Rev. date: 06/21